


Email : sales@equipmentsexporters.com Phone: +91-9311469084

	duct Code : RPHY10400002
--	-----------------------------

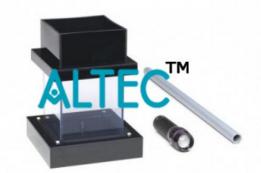

Description :

Cloud Chamber Cold Plate Type

Technical Specification :

This safe and affordable Cloud Chamber requires no dry ice or liquid nitrogen.

Instead, the apparatus uses a specially designed cold plate that is capable of retaining a low temperature for about one hour by freezing the plate in a standard freezer.

The cold plate can be used repeatedly by freezing it again.

The only consumables required to observe radiation tracks are ethanol (99.5%) and hot water (~50° C).

Tracks of radiation can be observed for about 20 minutes (until ethanol is vaporized) in a room at 25? (77°F).

Detailed Instruction Guide included.

Address: 75, Lajpat Nagar-IV, New Delhi-110024 Phone: +91-9311469084

